


THE ART AND SCIENCE OF DEFIBRILLATION

Lifeline VIEW

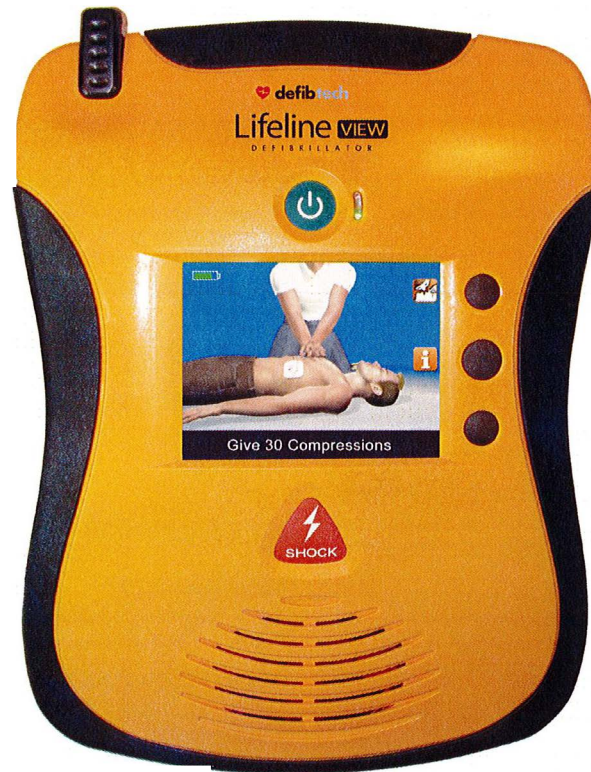

Lifeline VIEW Semi-Automatic Defibrillator

FULL-COLOR VIDEO INSTRUCTION DELIVERS BREAKTHROUGH EASE-OF-USE

The Lifeline VIEW makes it easy for anyone — from first responders to untrained bystanders — to respond confidently and appropriately during an emergency. This revolutionary life-saving device combines a sleek, user-friendly design with military-grade specifications — including the most proven biphasic waveform. It's also the first and only AED with a full-color interactive display that shows step-by-step videos for performing CPR, rescue breathing and external defibrillation.

You can also be confident that the device will be easy to maintain and ready for rescue. The Lifeline VIEW provides clear guidance for using and maintaining the device. Additionally, it has been designed to support upgrades and enhancements.

Like the original award-winning Defibtech AED, the Lifeline VIEW includes all the mission-critical features necessary to provide the most advanced cardiac treatment and exceeds the American Heart Association performance recommendations.



Bartec Fire Safety Systems Ltd.

8561 Government Street, Burnaby, B.C. V3N 4S9, Canada

604-420-2964 * sales@bartecfire.com * www.bartecfire.com

Lifeline VIEW Semi-Automatic Defibrillator

TECHNICAL SPECIFICATIONS†

DEFIBRILLATOR

TYPE

Semi-Automatic external defibrillator

MODEL

DDU-2300 series

WAVEFORM

Biphasic Truncated Exponential
(Impedance compensated)

ENERGY

Adult: 150 Joules
Child / Infant: 50 Joules
(Nominal into 50 Ohm load)

CONTROLS

Lighted On/Off button
Lighted Shock button

CHARGE TIME*

Less than 4 seconds
(shock decision to shock)

DISPLAY

High-resolution color LCD

VIDEO PROMPTS

Full motion video
On-screen text prompts

CPR COACHING

Video and voice coaching
On-demand video help

VOICE PROMPTS

Extensive voice prompts guide
user through operation of the unit

RESCUE PROTOCOL

AHA/ERC 2010
Supports protocol updates by
the user (password protected)

*typical, with new battery at 25° C

DEFIBRILLATION / MONITORING PADS

MODEL

Adult: DDP-2001
Child / Infant:
DDP-2002 and DDP-2003

SURFACE AREA**

Adult: 77 cm²
Child / Infant: 50 cm²

TYPE

Pre-connected, single-use,
non-polarized, disposable,
self-adhesive electrodes with
cable and connector

**nominal, each pad

EVENT DOCUMENTATION

INTERNAL EVENT RECORD

Critical ECG segments and rescue
event parameters are recorded
(greater than 60 minutes) and
can be downloaded to a removable
data card

PC-BASED EVENT REVIEW

ECG with event tag display, and
audio playback when available

REMOVABLE STORAGE

(optional) Up to 30 hours of ECG
and event data storage (no audio
option) or up to 3 hours of audio
(audio option). ECG and event
storage on a removable data card.
Actual length of storage is
dependent on card capacity.

USB PORT

Event download and
maintenance operations

PATIENT ANALYSIS SYSTEM

PATIENT ANALYSIS

Automatically evaluates patient
impedance for proper pad contact.
Monitors signal quality and analyzes
patient ECG for shockable/
non-shockable rhythms.

SENSITIVITY/SPECIFICITY

Meets AAMI-DF-80 specifications
and AHA recommendations

BATTERY PACK

MODEL

DBP-2003 (standard),
DBP-2013 (aviation; TSO C-142a)

POWER

12V, 2800 mAh

TYPE

Lithium/Manganese Dioxide
Disposable, recyclable,
non-rechargeable

CAPACITY*

125 shocks or 8 hours
continuous operation

STANDBY LIFE*

4 years

LOW BATTERY INDICATORS

Visible
Audible

*typical, with new battery at 25° C

ENVIRONMENTAL

TEMPERATURE

Operating: 0 to 50°C (32 to 122°F)
Standby: 0 to 50°C (32 to 122°F)

RELATIVE HUMIDITY

Operating / Standby: 5%-95%
(non-condensing)

ALTITUDE

-500 to 15,000 ft (-150 to 4500 m)
per MIL-STD 810F 500.4 Procedure II

VIBRATION

Ground (MIL-STD-810F 514.5
Category 20)

Helicopter (RTCA/DO-160D,
Section 8.8.2, Cat R, Zone 2,
Curve G)

Jet Aircraft (RTCA/DO-160D
Section 8, Cat H, Zone 2,
Curves B & R)

SHOCK / DROP ABUSE TOLERANCE

MIL-STD-810F 516.5 Procedure IV
48 inches (1.2 meters), any edge,
corner, or surface, in standby mode

CRUSH TEST

1,000 pounds (450 kg)

SEALING / WATER RESISTANCE

IEC 60529 class IP55;
Protected against dust and
waterjets (battery pack installed)

ESD

IEC 61000-4-2: (Open air up to
15kV or direct contact up to 8kV)

EMC (Emission)

CISPR 11 Group 1 Level B
and FCC Part 15

EMC (Immunity)

IEC 61000-4-3 and IEC 61000-4-8

SELF TESTS

AUTOMATIC

Automatic daily, weekly, monthly
and quarterly circuitry tests

BATTERY INSERTION

System integrity test on
battery insertion

PAD PRESENCE

Pads preconnected tested daily

USER-INITIATED

Unit and battery pack system
test initiated by the user

STATUS INDICATION

Visual and audible indication
of unit status

STATUS SCREEN

Unit self-test results
Pads and battery information
(status and expiration)

PHYSICAL

SIZE

7.3 x 9.5 x 2.3 Inches
(18.5 x 24 x 5.8 cm)

WEIGHT (with battery)

Less than 3 lbs (1.4 kg)



Bartec Fire Safety Systems Ltd.
8561 Government Street, Burnaby, B.C. V3N 4S9, Canada
604-420-2964 * sales@bartecfire.com * www.bartecfire.com



Defibtech, LLC • Guilford, CT 06437 USA • 1-203-453-4507 • 1-866-DEFIB-4U (1-866-333-4248)
www.defibtech.com

†Specifications subject to change without notice



The Lifeline VIEW AED
is cleared by the FDA

DAC-2703-BH

ELECTRONIC DISTRIBUTION